

Pick Up **New** Products

Product name: PFA bottle

Sanplatec Fluorocarbon Polymers Product Series(PFA bottle, Beaker, Volumetric Flask)



High airtightness without inner seal



The inner seal is made of the same PFA as the main body and cap.



Appearance inspection using a high-sensitivity camera

■ Features

- Bottles with superior heat resistance, cold resistance, chemical resistance, ultraviolet light resistance, non-stick, electrical insulation, and purity
- Temperature tolerance ranged from -196°C to +200°C.
- Extremely little eluate from the container, which does not contaminate the contents
- Strong for acids and alkalis and resistant to most organic solvents

■ Specification

- Leak (decompression) test and appearance test were performed.
- The total number of appearance examinations was performed.
- Controlled manufacturing lot
- Individual packaging
- Implement thorough quality control






■ Other Products

- Acid washed PFA bottle
- Clean packaged PFA bottle
- Pyrogen-free washed PFA bottle

Pick Up **New** Products

Product name: PFA bottle GL45, Volumetric Flask

Sanplatec Fluorocarbon Polymers Product Series(PFA bottle, Beaker, Volumetric Flask)

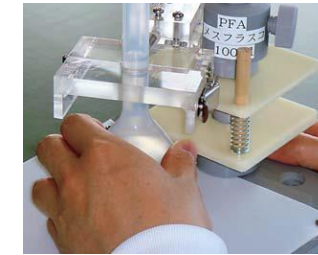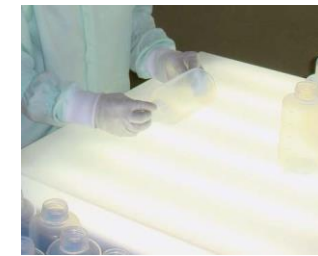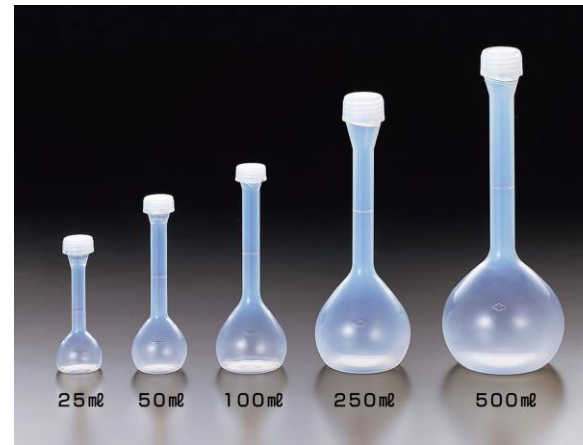

■ Features

- The mouth is compliant with GL45 and is compatible with other company GL45 caps. (Photo)
- The bottle body is thick
- Dimple Specification considering the liability of the 1L (Photo)

■ Features

- Resistant to most chemicals such as strong acids, strong bases and organic solvents.
- Since the elution of impurities is extremely small, it can be used for water quality inspection, reagent adjustment for ICP analysis, etc.
- Examination of appearance (total number of tests) was performed on an irradiated table to check for the presence of foreign matter.
- Weight is measured by manual operation of each piece, and the mark is stamped using the item cutting jig of the in-house developed product.

Pick Up **New** Products

Product name: PFA beaker, wash bottle, graduated cylinder

Sanplatec Fluorocarbon Polymers Product Series(PFA bottle, Beaker, Volumetric Flask)



■ PFA Wash Bottle ■



■ PFA Gas Wash Bottle ■



■ PFA Graduated Cylinder ■



■ Features

- The handle is also molded in one piece, so there is no fear of dislodging.
- Material marking mark, date of manufacture, etc. is stamped on the outer bottom of all sizes.
- Strong for acids and alkalis and resistant to most organic solvents